
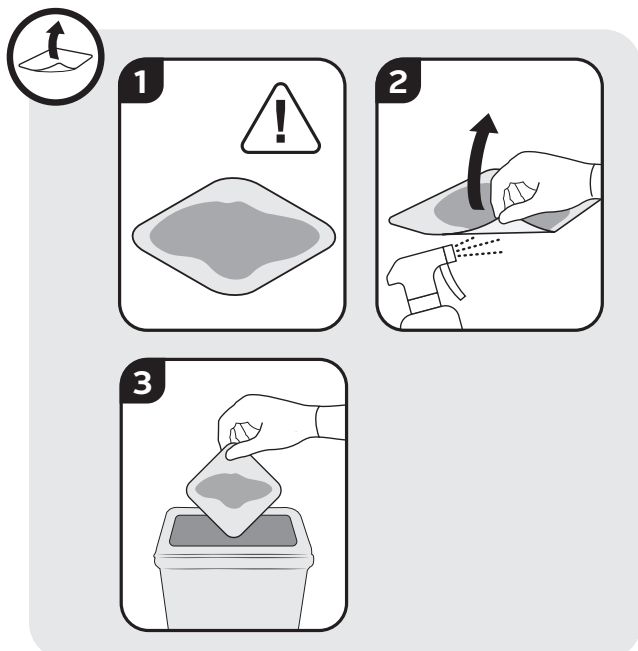
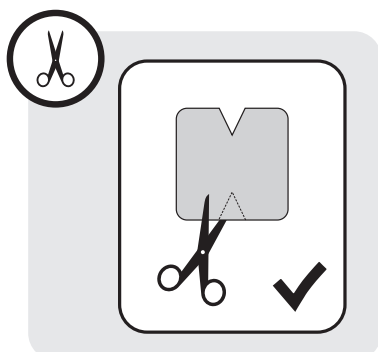
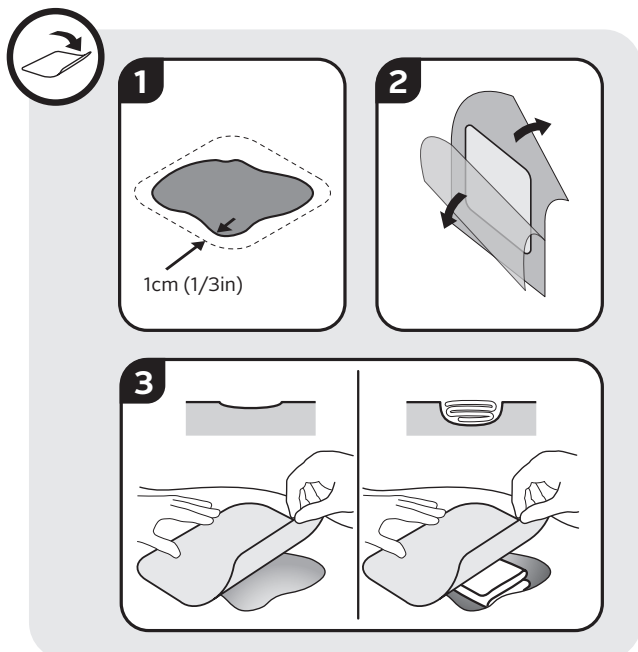
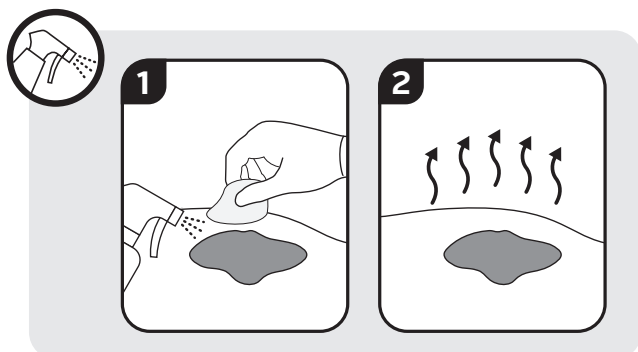
# ConvaFoam™

## Non-adhesive

### INSTRUCTIONS FOR USE / INSTRUCCIONES DE USO

Non Adhesive Hydrofiber™ Foam Dressing / Apósito no adhesivo de espuma y de Hydrofiber™
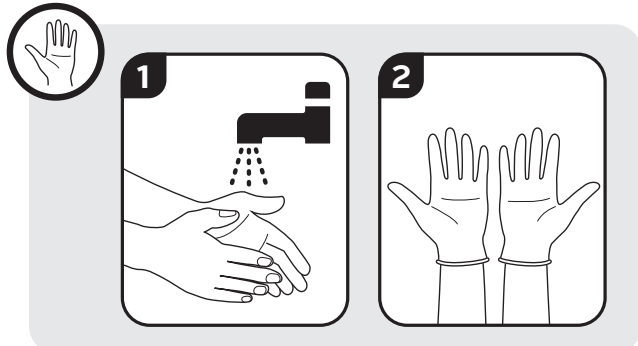
### PRODUCT DESCRIPTION

ConvaFoam™ Non-Adhesive is a sterile, Hydrofiber™ foam wound dressing, consisting of:

- a showerproof, breathable outer polyurethane film which provides a bacterial and viral barrier,
- a soft, multi-layered, absorbent central pad containing polyurethane foam, superabsorbent fibres and a Hydrofiber™ (Sodium carboxymethylcellulose) wound contact layer. The pad absorbs and retains moisture, exudate and bacteria.

### CLINICAL BENEFITS

ConvaFoam™ Non-Adhesive dressing is designed to provide a moist wound healing environment, aiding autolytic debridement and removing dead-space between the wound and dressing interface.

ConvaFoam™ Non-Adhesive dressing is designed to manage excess exudate levels which may further damage the wound bed and surrounding skin. ConvaFoam™ Non-Adhesive dressing is designed to protect intact skin from moisture, shear and friction damage.

### INTENDED USE

ConvaFoam™ Non-Adhesive dressings may be used, with the consultation of a healthcare professional, for the management of exuding and non-exuding wounds, and for the protection of intact skin.

ConvaFoam™ Non-Adhesive dressings may also be used on intact skin as part of a pressure injury prevention protocol.

ConvaFoam™ Non-Adhesive dressings may be used as a primary or secondary dressing.

### INDICATIONS

ConvaFoam™ Non-Adhesive dressing is indicated for the management of:

- Leg Ulcers
- Pressure Injuries
- Diabetic Foot Ulcers
- Surgical Wounds
- Partial Thickness Burns
- Traumatic Wounds

### CONTRAINDICATIONS

ConvaFoam™ Non-Adhesive dressings should not be used on individuals who are sensitive to or who have had an allergic reaction to the dressing or components stated in the Product Description above.

### INTENDED USERS

ConvaFoam™ Non-Adhesive dressings are used by health care professionals, patients and caregivers.

### TARGET PATIENT POPULATION

ConvaFoam™ Non-Adhesive dressings are intended to be used on patients with one of the wound types listed in the indications, those who require skin protection from breakdown or for the prevention of pressure injury on intact skin as part of a pressure injury prevention protocol.

### PRECAUTIONS

ConvaFoam™ Non-Adhesive dressings are for single-use only. Do not reuse.

Reuse may lead to increased risk of infection, cross-contamination and delayed healing.

ConvaFoam™ Non-Adhesive dressings are not compatible with oil-based products.

Do not apply oxidising solutions such as hydrogen peroxide or hypochlorite solutions to the dressing.

If you are using oxidising agents as part of a protocol of care, ensure that the wound and surrounding skin is dry and free from any residues of oxidising agents before applying the dressing.

For highly exuding wounds, ConvaFoam™ Non-Adhesive dressings may be used as a secondary dressing in conjunction with a primary dressing, such as Aquacel™ silver and non-silver Hydrofiber™ products.

ConvaFoam™ Non-Adhesive dressings should not be submerged during bathing or swimming.

Consult a healthcare professional if any of the following are observed:

- Irritation (reddening, inflammation)
- Maceration (whitening of skin)
- Hyper-granulation (excess tissue formation)
- Sensitivity or allergic reaction
- Signs of infection (increased pain, bleeding, warmth/redness of surrounding tissue, wound exudate)
- Changes in wound colour and/or odour
- Increase in wound size
- Wound not showing signs of healing.

After use, this product may be a potential biohazard. Handle and dispose of in accordance with accepted medical practice and applicable local laws and regulations.

If during the use of this device, or as a result of its use a serious incident has occurred, please report it to the manufacturer at [www.convatec.com](http://www.convatec.com) and to your national authority.

### CONTACT AND DURATION OF USE

ConvaFoam™ Non-Adhesive dressings can be worn for up to 7 days, dressings should be changed earlier if clinically indicated, and are safe for continuous use.

### DIRECTIONS FOR USE

#### 1. Storage

- Store at room temperature.
- Protect from light.
- Keep dry.
- Ensure the dressing pouch is not damaged or opened prior to use to guarantee sterility. If the dressing or pouch is damaged or opened, do not use and dispose according to local regulations.
- Do not use after expiry date.

#### 2. Wound Site Preparation and Skin Cleansing:

- Before applying the dressing ensure the skin

surrounding the wound is free from creams, moisturisers and oil-based products.

- Before applying the dressing cleanse with a suitable wound cleanser and dry the surrounding skin.
- If a moisture barrier is required, select a spray or liquid that will dry before dressing application.

### 3. Dressing Preparation and Application:

- When applying to a wound choose an appropriate dressing size where the central absorbent pad is at least 1cm (1/3 inch) larger than the wound.
- When applying to intact skin choose a dressing which covers the entire vulnerable area.
- Remove the dressing from the sterile pack.
- Hold the dressing over the wound or area of skin to be protected and line up the centre of the dressing with the centre of the wound or area of skin to be protected.
- Ensure the dressing is not stretched upon application. Ensure joints are sufficiently flexed to allow patient mobility.

If required, dressings may be cut to fit. Discard any unused portion of the product after dressing the wound. Use additional tape or film dressing to create a seal.

### 4. Dressing Removal:

- The dressing should be changed when clinically indicated, when saturated with wound fluid, if it leaks, or if it becomes damaged or soiled.
- If removing the dressing is difficult, the dressing should be soaked with sterile saline or water and removed slowly.
- To remove the dressing, press down gently on the skin immediately surrounding the dressing and carefully lift one corner of the dressing. Continue lifting the dressing until all edges are free. Carefully lift away the dressing and discard as per local clinical protocols.

### Protective Care:

During protective care only, the dressing corner may be lifted for inspection of skin and resealed. The dressing must be changed if skin is broken.

If further information or guidance is needed, or to report an adverse event, please contact Convatec, contact details can be found at [www.convatec.com](http://www.convatec.com) or use the contact details at the end of this document.

Made in UK

© 2024 Convatec Inc.

™ indicates a trademark of Convatec Inc.

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Convatec Limited